

## Path-Chek® Hygiene protein

A rapid test for the detection of food residues from environmental surfaces or processing equipment.

- Results in 5 seconds
- No Laboratory facilities required
- Simple to use
- Cost effective - 2 tests per swab

Path-Chek® Hygiene Protein is an inexpensive and rapid chemical based test for the detection of food residues on food manufacturing surfaces after routine cleaning and sanitising procedures have been performed. Food manufacturing surfaces must be free of any food residue to be considered clean, as the efficient removal of food residues by applying appropriate cleaning and sanitising procedures will also remove any bacteria present on the surface.



Path-Chek® Hygiene Protein is based on the detection of invisible food residues remaining on surfaces when cleaning and sanitising procedures are not effective or have not been performed correctly. The test provides an instant result without the need for laboratory facilities or incubation. Path-Chek® Hygiene Protein does not detect bacteria so tests for specific food borne pathogens should always be performed but protein and food residues detected by the test indicate that bacteria may also be present on the surface.


Path-Chek® Hygiene Protein swabs are economical to use as one swab has two swab heads and therefore two tests per swab. If you only require to use one test then you can store opened swabs in a zip lock bag for up to 7 days at 2-8°C.

Easy to use  
Cost effective

# Path-Chek<sup>®</sup> Hygiene protein

How to use Path-Chek<sup>®</sup> Hygiene protein:

1. Separate the number of swabs required.
2. Tear the swab at the mid-point as indicated.
3. Remove one side of the foil package to expose a swab head.
4. Test sample area with exposed swab head.
5. A positive reaction is indicated by the swab head changing from yellow to green within 5 seconds.
6. Use removed foil package to cover used swab head.

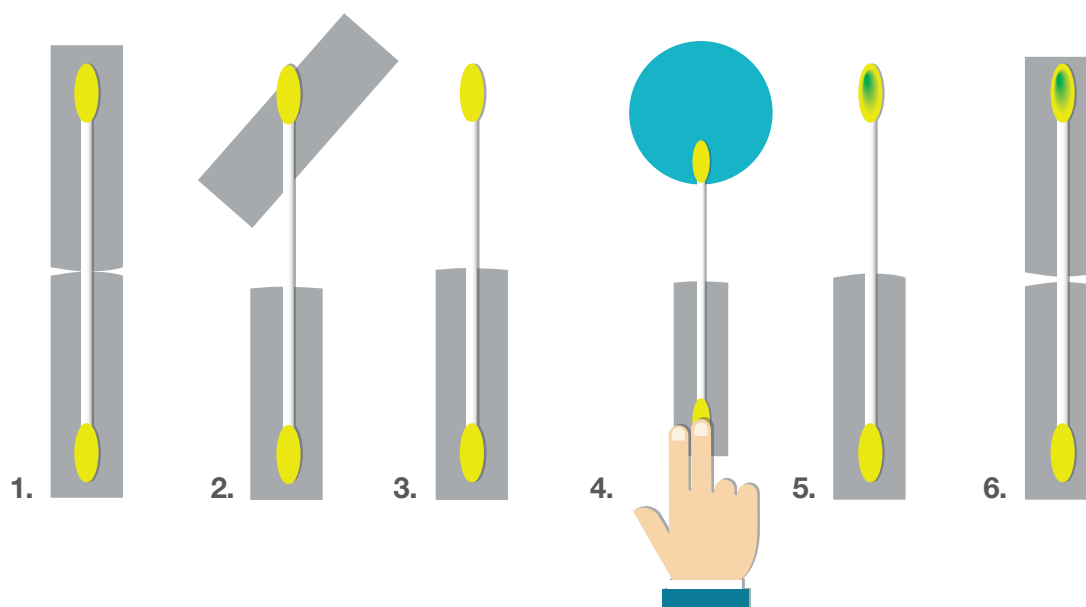

## Ordering Information:

Code	Description	Size	Storage
PC 006	Path-Chek <sup>®</sup> protein	50 swabs (100 tests)	Room temperature

If you need any further information on any of the above tests please contact Microgen Bioproducts Ltd.

**MICROGEN**  
BIOPRODUCTS

**Microgen Bioproducts Ltd**, 1 Admiralty Way, Camberley, Surrey GU15 3DT, United Kingdom  
T: +44 (0)1276 600 081 F: +44 (0)1276 600 151 E: sales.camberley@lab21.com  
[www.microgenbioproducts.com](http://www.microgenbioproducts.com)

GMM-106-EN revision 1.

Part of the  
**NOVACYT**  
GROUP